

Tierrett Small Foley Catheter

Indwelling Bladder Catheter

Registration No. of Medical Equipment for Animal Use: 28 Douyaku 2048

Developed for use in **“very small dogs”**

“4.5Fr” silicon balloon catheter has been made available through a combination of various technologies.



○ Features

- Very small-size and can be smoothly inserted even in very small dogs
- Comes with hydrophilic Ni-Ti stylet
- Comes with multi-cap that can be used in 3 ways


“Tierrett”, derived from the German word **“Tierrettung”** (aid for an animal), is a brand of **CREATE MEDIC** for animal medical products that were developed with the aim of helping sick animals.



Tierrett Small Foley Catheter

Feature 1.

- Achieves both insertability and removability
Comes with moderately elastic hydrophilic stylet (with Ni-Ti cable core)
 - High insertability is ensured by a hyper-elastic Ni-Ti (nickel titanium alloy) cable core.
 - Removability of the stylet is also ensured by an outer hydrophilic coating.

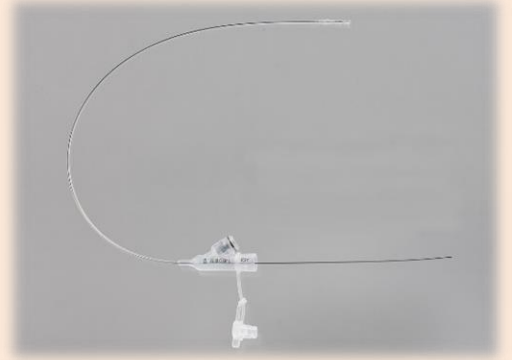

Feature 2.

- Useful **multi-cap** has both syringe connection and plugging functions.
 - It can be used in 3 ways as a catheter tip, lure tip, and plug.



Feature 3.

- Easy-to-handle **lightweight/compact joint**
 - Reduced weight to 50% of conventional company products and 55% of length

Feature 4

- Made of human-body-friendly **silicon**
 - Hypoallergenic to urethral mucosa, reducing the adherence of urine constituent

Product Specification

Size designation	Outer diameter	Balloon capacity	Effective length	Product code	JAN code
4.5Fr	1.5mm	0.8mL	300mm	800-001-4702	4542187203738
			550mm	800-001-4703	4542187203745
5Fr	1.7mm	0.8mL	300mm	800-001-4704	4542187203752
			550mm	800-001-4705	4542187203769

5-pack box/EOG sterilized

- Please note that the product specifications are subject to change without notice to reflect the results of our consistent R&D efforts.
- Read the package insert before using the product.



CREATE MEDIC CO.,LTD.

<https://www.createmedic.co.jp/english/>

【Contact information】

Animal Medical Business Department E-mail: oem-infoch@createmedic.co.jp

【Head Office】 2-5-25 Chigasaki Minami, Tsuzuki-ku, Yokohama, Kanagawa, Japan

